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### Instruction Manual

---

Thank you for choosing our product. We hope you will enjoy using the appliance.

#### Symbols in these Instructions for Use

Important information for your safety is specially marked. It is essential to comply with these instructions in order to avoid accidents and prevent damage to the machine:

#### **WARNING:**

This warns you of dangers to your health and indicates possible injury risks.

#### **CAUTION:**

This refers to possible hazards to the machine or other objects.

#### **NOTE:**

This highlights tips and information.

### General Safety Instructions

---

Read the operating instructions carefully before putting the appliance into operation and keep the instructions including the warranty, the receipt and, if possible, the box with the internal packing. If you give this device to other people, please also pass on the operating instructions.

- The appliance is designed exclusively for private use and for the envisaged purpose. This appliance is not fit for commercial use.

- Do not use it outdoors. Keep it away from sources of heat, direct sunlight, humidity (never dip it into any liquid) and sharp edges. Do not use the appliance with wet hands. If the appliance is humid or wet, unplug it immediately.
- When cleaning or putting it away, switch off the appliance and always pull out the plug from the socket (pull the plug itself, not the lead) if the appliance is not being used and remove the attached accessories.
- Do not operate the machine without supervision. If you leave the room you should always turn the device off. Remove the plug from the socket.
- The device and the mains lead have to be checked regularly for signs of damage. If damage is found the device must not be used.
- Do not try to repair the appliance on your own. Always contact an authorized technician. To avoid the exposure to danger, always have faulty cable be replaced only by the manufacturer, by our customer service or by qualified person and with a cable of the same type.
- Use only original spare parts.
- In order to ensure your children's safety, please keep all packaging (plastic bags, boxes, polystyrene etc.) out of their reach.

#### **WARNING!**

**Caution!** Do not allow small children to play with the foil as there is a **danger of suffocation!**

- Pay careful attention to the following "Special Safety Instructions".

## Special safety instructions for this device

---

#### **WARNING:**

- This device can be used by children over the age of 8 and by persons with reduced physical, sensory or mental capabilities or lack of experience and/or knowledge if they are supervised.
- **Children** may not play with the device.
- The device may only be used with cold or lukewarm water.
- Only use the device for the cleaning of the mouth.
- Empty the water container after each application to prevent an accumulation of bacteria.

- When using braces, find out whether the usage of an oral irrigator is possible. We do not assume any liability for possible damages.
- In case of diseases in the mouth and throat area, you should ask your doctor whether the use of an oral irrigator is advisable.
- Ensure that the water tube (connection of device and oral irrigator) is not stuck or kinked.

## Overview of the control elements

- 1 Single jet nozzle
- 2 Water stop switch
- 3 Handle
- 4 Water hose
- 5 Water pressure regulator
- 6 On/Off switch
- 7 Unit base
- 8 Water container
- 9 Handle for opening the lid / Handle for water container

Level I:	low water pressure
Level II:	medium water pressure
Level III:	high water pressure

### Water stop switch (2)

You can use this switch to momentarily interrupt the application. The shut-off is done mechanically. The pump of the device will continue to be in operation.

Position I of the switch:	Water is flowing
Position 0 of the switch:	Water is stopping

### Nozzles (1)

4 nozzles with different color markings are contained in the scope of delivery. These can be assigned to different users.

### When the gums are bleeding . . .

It is normal for a normal bleeding to occur when using the oral irrigator the first few times. Should your gums bleed continually, please consult your dentist.

## Unpacking the Device

1. Take the device out of its packaging.
2. Remove all packaging materials, such as foils, filler material, cable brackets and cardboard packaging.
3. Check the scope of delivery.

### NOTE:

There may still be some dust or production residues on the device. We recommend that you clean the device as described under "Cleaning".

## Application instructions

### Electrical connection

1. Check whether the mains voltage that you want to use matches the one of the device. You will find the information about it on the type plate.
2. Connect the mains cable to a properly installed grounded socket.

### On / Off switch (6)

Switch setting 0:	The device is turned off
Switch setting I:	The device is turned on

### Water pressure regulator (5)

With the water pressure regulator, you can set three different power levels.

## Initial operation

1. Remove the water container (8) from the unit base by the handle. Open it with the lever (9).
2. Fill the water container with water. Observe the marking MAX. Close the water container. Set it back onto the unit base.
3. Turn on the device with the On / Off switch (6).
4. Remove the handle (3) of the unit base. Hold the opening above a wash basin.
5. Switch the regulator (5) to level III.
6. Unlock the water stop switch (2) by sliding it to position I.
7. Let the water run directly into the wash basin to remove possible residues. Then lock the water stop switch by sliding it to the position 0.
8. Turn off the device.
9. Briefly rinse off the nozzles under running water.

Your oral irrigator is now ready for use.

## 24 English

### Usage

Before using the water jet you should clean your teeth with a toothbrush.

1. Insert the nozzle onto the handle as far as possible. Turn the nozzle so that its opening points to the water stop switch. Ensure that the nozzle is fixed firmly in position.
2. Fill the water container with water which is hand warm. Do not fill it with hot water!
3. In order to clean the nozzles, allow water to flow through the nozzle head (at maximum water pressure) into the wash basin before and after each time the device is used.
4. Then use the water pressure regulator to select the desired water pressure. Begin on the smallest setting, especially if you have problems with bleeding of the gums or when the water jet is used by children.
5. Unlock the water stop switch on the handle.
6. Lean over the wash basin. Guide the water jet along the gums, around the teeth.
7. Lock the water stop switch after the application. Set the water pressure regulator to level I.
8. Turn off the device.
9. Please clean the device as described under "Cleaning and Care".

### Cleaning and Care

#### WARNING:

- Pull the mains plug from the socket before each cleaning.
- Do not immerse the unit base and the mains plug in water.

#### CAUTION:

- Do not use any wire brush or other abrasive objects.
  - Do not use any sharp or abrasive cleaners.
1. Pour out the remaining water from the water container.
  2. Pull the nozzle off the handle. Clean the nozzles under warm running water. Dry the nozzles and place them onto the unit base.
  3. If you have used water additives, rinse the water container with clear water. Then let ca. 300 ml clear water run through the device as described under "Initial operation".
  4. Dry the handle before you place it back onto the unit base.
  5. Dry off the water container.
  6. Clean the unit base with a moist towel. Then wipe it dry.

### Noise development

The emission value referred to the place of work is below 70 dB (A).

#### Operating noises


The noise that you perceive during use occur by the water flow. These noises are unavoidable due to the high performance, but make no shortage represents.

### Technical Data

Model: ..... MD 5613  
Power supply:.....230 V~ ; 50 Hz  
Power consumption:.....max. 46 W  
Protection class:.....II (IPX 7)  
Net weight:..... 1.3 kg

The right to make technical and design modifications in the course of continuous product development remains reserved.

This device has been tested according to all relevant current CE guidelines, such as electromagnetic compatibility and low voltage directives, and has been constructed in accordance with the latest safety regulations.
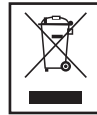
### Disposal

#### Meaning of the "Dustbin" Symbol

Protect our environment: do not dispose of electrical equipment in the domestic waste.

Please return any electrical equipment that you will no longer use to the collection points provided for their disposal.

This helps avoid the potential effects of incorrect disposal on the environment and human health.

This will contribute to the recycling and other forms of reutilization of electrical and electronic equipment.

Information concerning where the equipment can be disposed of can be obtained from your local authority.

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