
## Instruction Manual BS-8050

### Bathroom scales
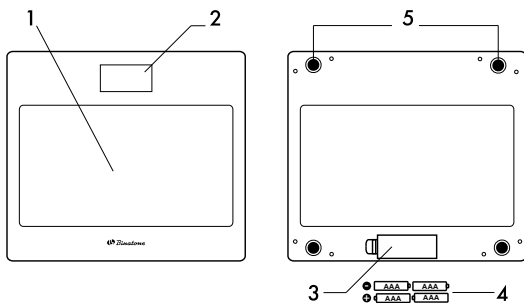
## Bathroom scales BS-8050

The appliance is designed for static human body measurement at home.

### Technical specifications and features:

Max. weight: 150 kg

Deviation: 0.1 kg

Power source: 4 batteries 1.5V type AAA

### SAFETY MEASURES

This appliance is powered by batteries, therefore during operation it is necessary to comply with the following safety measures:

- ▶ If the appliance is not going to be used for a long period of time, it is recommended to take the batteries out of the appliance.
- ▶ Only use the batteries as recommended by the manufacturer.

### SPECIAL SAFETY MEASURES

- ▶ Only perform the weighing of children and people with limited physical abilities under strict supervision.
- ▶ Do not stand on the platform with wet feet since it can be slippery. Always make sure that the scales are dry.
- ▶ While weighing, always stand in the central part of the platform in order not to overturn the scales or fall off.
- ▶ Store the appliance in the original package, in a horizontal position.

### DESCRIPTION OF PARTS OF THE APPLIANCE

1. Platform
2. LCD display with blue backlight
3. Battery compartment
4. 4 batteries 1.5V type AAA
5. Rubber feet

### BEFORE THE FIRST USE

Open the battery compartment (3) on the platform (1) and remove the insulation to release contacts.

### OPERATION

1. Put the scales on a flat solid horizontal surface. Do not put the scales on carpets or other soft covering, since this will affect the measurement accuracy.
2. To turn on the scales, stand onto the scales with your foot evenly distributing your weight over the platform (1) and don't move. Wait a few seconds until the appliance calculates the result, which will appear on the display.
3. If the measured weight exceeds the allowable load (150 kg), the display will show the indication of excess load «O-Ld».
4. If «LO» appears on the scales display (2), this means the battery requires replacement. Put the new battery in the appliance, observing polarity.

### CARE AND CLEANING

1. Wipe the appliance with a damp lint free cloth, wiping it dry, so that there are no water marks on the platform (1). Do prevent water getting inside the appliance. Do not use hard abrasive materials and substances for cleaning.
2. Always keep the appliance in a horizontal position in a dry, cool place. Do not store the appliance in places with high humidity or temperature.

### AFTER-SALES SERVICE

If your appliance does not operate or does not operate properly, contact an authorized service center for consultancy or repair.

*The warranty does not cover:*

- product packing, consumables and accessories (knives, tips, etc.);
- mechanical damages and natural wear of the product;
- failures due to overloading, improper or careless maintenance, penetration of liquids, dust, insects and other foreign matters inside the product, exposure of plastic and other not thermal-resistant components to high temperatures, force-majeure circumstances (accidents, fire, flooding, power grid failure, lightning stroke, etc.);
- damages caused by:
  - ▶ Non-fulfillment of the above-mentioned precautions and operating manuals.
  - ▶ Application in commercial purposes (application of the product beyond personal daily living needs).
  - ▶ Unauthorized repair or replacement of the components of the appliance, performed by an unauthorized person or service center.

Also the warranty doesn't cover replaceable batteries.

Manufacturing date is specified on the appliance.

RU XX YY ZZZ code reads as:

XX – production year

YY – production month

ZZZ – batch number

**Major Importer of Binatone TM in Russia:** LLC TORGPROEKT Domodedovo, Address of the Company: house 7/1, Pogorelskiy per, Moscow, 119017

**Manufacturer:** Binatone Industries Ltd, Great Britain, U.K., 10th Floor, 6 Mitre Passage, Greenwich, London, SE 10 0ER.

Made in PRC for Binatone Industries Ltd, Great Britain.

**Warranty** 2 years.

**Exploitation period:** 3 years after the day of purchase.

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