Elenberg

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CONTENTS

THE DEAR BUYER!

Congratulations, you made a good choice!

You became the owner of ELENBERG production, which offers progressive design and high quality. For full and safely usage of the given device, please, closely read the instruction manual.

USER-FRIENDLY TIPS

During Measuring

- A. Always remove your shoes and socks, wipe clean your feet before using.
- B. It is advisable to take measurement at the same time of the day.
- C. Measured result could be misleading after intensive exercise, excessive dieting or under extreme dehydration condition.
- D. Always weigh/use the scale on a hard and flat surface. When you are measuring, please keep still.
- E. The body fat/hydration measured data of the following people may have deviation:
 - Children under 10 or adult over 100 (can only use the device under normal weighing mode)
 - Adult over 70
 - Body building athletes or other occupational athletes.
- F. People with the following conditions should not use the device:
 - Pregnant women
 - Those on dialysis treatment
 - Those who use heart pacemaker or with other implanted medical device.

Safety Warning

- A. The Fat% Hydration% Muscle% content evaluated by this device is only for reference (not for medical purpose). If your Fat% Hydration% content is over/below normal level, please consult your doctor for more advice.
- B. For the data accuracy of the scale, please keep it away from strong electromagnetic field.
- C. It could be slippery if you weigh on a wet surface. Moreover, if the weighing platform surface is glazed, never step on the edge of the scale and do keep balance during your weighing.

FEATURES

- 8 mm tempered glass platform
- Ultra slim design: thickness 20 mm
- Measure body content of: % fat, % hydration and % muscle
- Memory for 12 persons
- Max weight: 150 kg
- Division: 100 g
- Touch button
- Auto on / off
- Power: 2 x 3 V CR2032 lithium cells (included)


PACKING LIST

- Scale
- Two batteries "CR2032"
- Instruction manual
- Warranty card

BATTERY

Please, use 2 x 3 V CR2032 lithium batteries. Remove isolation sheet before use; Open the battery case cover at the bottom of scale, and install or replace battery as follow:

- Take out the used battery with the aid of a sharp object, as shown in picture.
- Install new battery by putting one side of the battery beneath the battery baffle and then pressing down the other side.


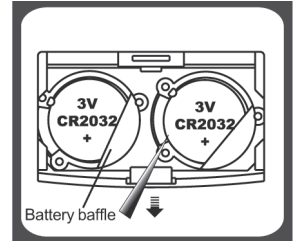
TO SET WEIGHT UNIT

If your scale is available (or applicable) in kg/lb/st, you can set to the one you prefer by pressing the unit conversion button (picture on the right). The unit conversion button locates on the bottom of scale.

Press the unit conversion button when scale is power off, LCD shows current weight unit.


Press the button again to convert weight unit in following sequence:

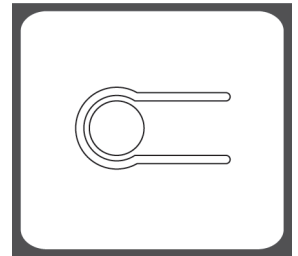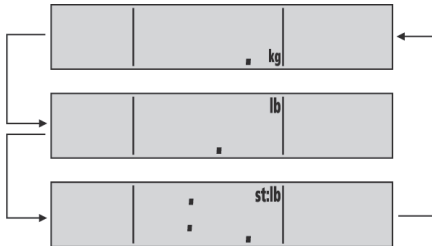After selecting the preferred unit, the scale will be auto-off after 3 seconds. Selected weight unit will be activated when the scale is turned on again.

- When setting lb/st as weight unit, the default height unit is feet/inch;
- When setting kg as weight unit, the default height unit is cm.

LCD DISPLAY

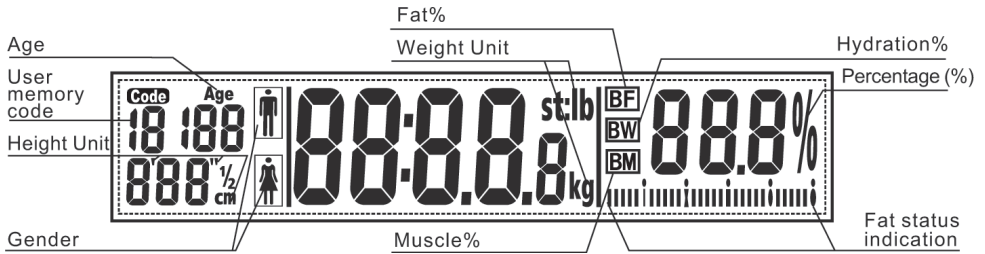KEY ILLUSTRATION

- set** Turn on the scale and prompt to next parameter. (Please refer to the Setting Personal Parameter over-leaf)
- ▲ Advance the value by one interval each press. Continuous pressing will speed up the advance of value.
- ▼ Decrease the value by one interval each press. Continuous pressing will speed up the decrease of value.

NORMAL WEIGHING MODE

Auto step on weighing


1. Install battery, and place the scale on flat floor.
2. Step on the weighing platform.
3. Digits flash, stabilized and locked, indicating your weight.

Consecutive add-on weighing


After digits stabilize, more additional loadings can be conducted. Ensure each add-on loading has to exceed 2 kg to get a new reading.

Example: If you want to get a weighing value of your baby, you can weigh yourself first, and then hold your baby together, without stepping off the scale. Then you will get the total weight reading. Subtracting the previous reading of yourself from the new reading, you will get the difference as your baby's weighing value.

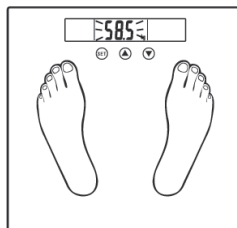Digits flash, stabilized and locked

Hold your baby without stepping off the scale.

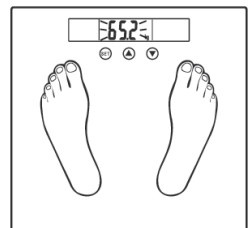Display total weight.

FAT/HYDRATION/MUSCLE MEASURING MODE


To set personal parameters

- At the first time when you install battery and switch on, all default parameters are set as follow: Gender - Male; Height - 165 cm; Age - 25.
- Parameters can be set in the following range:

Memory	Gender	Height	Age
1~12	male/female	100~250 cm	10 ~ 100

- Parameter setting

Example: the 5th user, female, 166 cm, 26 years old can set her parameters as follow:

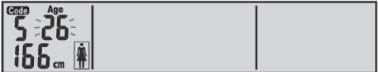A. Press "SET" to switch on the unit.	
B. Press ▲/▼ to select user code.	
C. Press SET to confirm and switch to next setting.	
D. Press ▲/▼ to select gender.	
E. Press SET to confirm and switch to next setting.	
F. Press ▲/▼ to select height.	
G. Press SET to confirm and switch to next setting.	
H. Press ▲/▼ to select age. Slightly tap on the scale platform by your hand or your foot. (Your personal parameters have been saved).	
I. LCD shows last measuring record of the 5th user (No record in the first time of use).	

Note:

- During the parameter-setting, you can always tap on the weighing platform any time to set the scale to zero-point.
- You can start measuring body fat/hydration based on the parameters just set/updated into the scale.
- After completing the setting cycle, you can always continue on pressing again to go through another setting cycle of all parameters.

To Start Measuring


To ensure accuracy, please remove your shoes & socks before stepping on scale.

- You can set Your parameters, as described in previous chapter or load al-

ready saved parameters. To load Your parameters press "SET" button and select user code by pressing ▲/▼ buttons.

Example: the 5th user, female, 166 cm, 26 years.

2. Step on scale.
4. Wait until the digits stabilize.
Weight readings locked and flash for 1 time.
5. After the weight is fixed, scale will begin to measuring Fat/Hydration/Muscle.

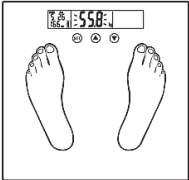6. The result will be displayed 3 times in following sequence:

	Fat % and status
	Hydration %
	Muscle %

The "||||| ||||| |||||" indication means that You are in good shape.

WARNING INDICATIONS

Low Battery Indication The battery power is running low, please replace with a new battery.	
Over-load Indication The weighing subject on the platform exceeds the maximum capability of scale. Please step off to avoid damage.	
Low Fat% Indication The Fat% is too low. Please lead a more nutritional diet and take good care of yourself.	
High Fat% Indication The Fat% is too high. Please watch your diet and do more exercise.	
Re-test Error exists, please re-test to get the correct result.	

REFERENCE INFORMATION

About Fat%

Fat% is an index % of fat content in human body.

The fat content in human body is of a certain level, overabundance of fat will bring bad effects on metabolism, thereby, various illnesses may invade.

Here goes an old saying, Obesity is hotbed of illness .

Working Principles

As electrical signal of certain frequency transmits through human body, the impedance of fat is stronger than that of muscle or other human tissues.

Through the use of biological resistance analysis, this product will pass a safe frequency of low voltage electrical current through human body, thus, the Fat% will be reflected according to the resistance differences when the electrical signal transmitting through human body.

Fat/Hydration Content Fitness Assessment Chart

After your weighing, you can use the following Assessment Chart for reference:

Age	Fat (Female)	Hydration (Female)	Fat (Male)	Hydration (Male)	Status	Key Line Indication
≤30	4.0 - 16.0%	66.0 - 57.8%	4.0 - 11.0%	66.0 - 61.2%	Too lean	
	16.1 - 20.5%	57.7 - 54.7%	11.1 - 15.5%	61.1 - 58.1%	Lean	
	20.6 - 25.0%	54.6 - 51.6%	15.6 - 20.0%	58.0 - 55.0%	Normal	
	25.1 - 30.5%	51.5 - 47.8%	20.1 - 24.5%	54.9 - 51.9%	Fat	
	30.6 - 60.0%	47.7 - 27.5%	24.6 - 60.0%	51.8 - 27.5%	Too fat	
>30	4.0 - 20.0%	66.0 - 55.0%	4.0 - 15.0%	66.0 - 58.4%	Too lean	
	20.1 - 25.0%	54.9 - 51.6%	15.1 - 19.5%	58.3 - 55.3%	Lean	
	25.1 - 30.0%	51.5 - 48.1%	19.6 - 24.0%	55.2 - 52.3%	Normal	
	30.1 - 35.0%	48.0 - 44.7%	24.1 - 28.5%	52.2 - 49.2%	Fat	
	35.1 - 60.0%	44.6 - 27.5%	28.6 - 60.0%	49.1 - 27.5%	Too fat	

* - The above information is only for reference.

CLEANING & MAINTENANCE

1. Use a slightly damp soft fabric to wipe clean the scale. DO NOT immerse it into water. Any use of chemical or corrosive detergent is not acceptable.
2. If the plastic components of the scale contacts with grease, spicy sauce, vinegar or other irritating substance, please clean it at once (Please keep the scale away from any fruit juice).
3. Always weight/use the scale on a hard and flat surface. DO NOT use on carpets or mats.

PERFORMANCE INDEX

Max Weighing Capacity	150 kg
Division	0.1 kg
Fat% Range	4.0% ~ 60.0%
Fat% Division	0.1 %
Hydration% Range	27.5% ~ 66.0%
Hydration% Division	0.1 %
Muscle% Range	20.0% ~ 56.0%
Muscle Division	0.1 %

DESIGN AND SPECIFICATIONS ARE SUBJECT TO CHANGE WITHOUT NOTICE!

Address to service center for consultation, when the service life of product has expired. Otherwise the further operation can entail impossibility of normal use of the product.

Service life of the given product - 3 years from the moment of sale