BS1800

WAGA OSOBOWA FITNESS

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Congratulations on your choice of our appliance, and welcome among the users of Zelmer products.

Please read this instruction manual carefully. Pay special attention to important safety instructions. Keep this instruction manual for future reference.

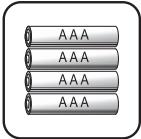
Handling instructions

Do not hit

Do not use or store in humid locations

Remove batteries if appliance is not to be used for a long period of time

Avoid overloading

Do not place any objects on the surface of the appliance

Do not immerse in water

Be careful not to slip when you step on the scale

People with pacemakers should not use the scale. In case of doubts ask your doctor for advice.

Pregnant women should use the weight reading function only.

The appliance is intended solely for household use, indoors.

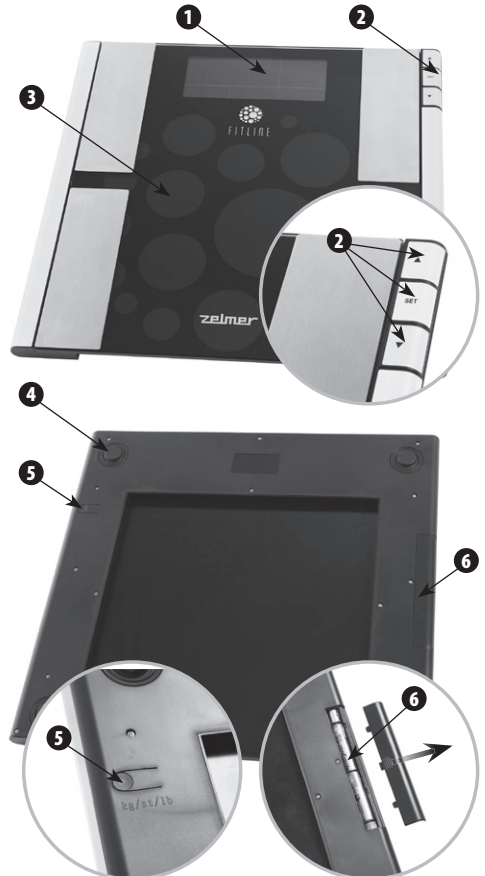
Technical data

Technical parameters are stated on the rating plate. ZELMER scales meet the requirements of relevant standards. The product is CE marked on the rating plate.

- Maximum load: 150 kg / 330 lb / 23 st 13 lb.
- Graduation: 0,1 kg / 0,2 lb.
- Battery type: 1,5 V AAA.


Description of the appliance

- 1 Electronic display
- 2 Function buttons
- 3 Weighing platform with BIA feature (bioelectrical impedance analysis)
- 4 Rubber feet
- 5 Weighing unit button
- 6 Battery compartment (4 x 1.5 V AAA)

Control panel

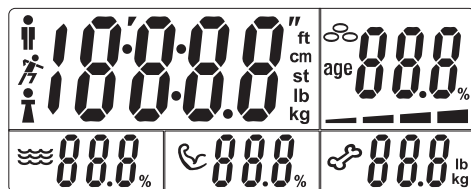DISPLAY SYMBOLS

- gender and body type
- body water
- muscle mass
- bone mass
- body fat
- age

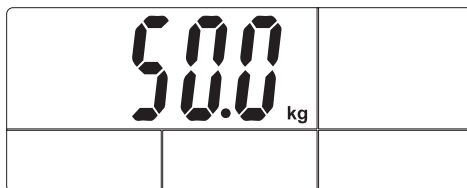
FUNCTION BUTTONS

- SET – set / confirm
- ▲ – increase / select
- ▼ – decrease / select

Use and functions

HOW TO USE

- ① Remove the scales from the box.
- ② Insert batteries into the battery compartment (see "REPLACING BATTERIES").
- ③ Switch on the scales by stepping on weighing platform.
- ④ Press the weighing unit button repeatedly until your required weighing unit shows on the display. The weighing unit button is located on the underside of the appliance.
- ⑤ Enter personal data of the user (see „BIA FEATURE“).
- ⑥ To start your measurement, wait until the display shows **00.0**. Step on the weighing platform with feet on the metal plates.
- ⑦a If the reading has stabilised but the user has not been identified, the result will be displayed for 5 seconds and the appliance will switch off automatically (weight measurement function).

i The user will not be identified if feet are isolated from the sensors on the weighing platform (e.g. by shoes).

- ⑦b The measurement will start once the reading has stabilised, the personal data of the user has been entered and the user has been identified.

In addition to displaying the body weight, the device provides information on the percentage of bone, muscle, water and fat in the body.

i This information will appear on the display for about 20 seconds and then the device will switch off automatically.

Compare your results with the tables below (see **Reference tables**).

REPLACING BATTERIES

The battery compartment is located on the underside of the appliance.

- ① To replace the batteries, pull out the compartment cover and remove old batteries.
- ② Insert 4 new AAA 1.5 V batteries.

! Make sure to insert the batteries the right way round.

i Check the inside of the compartment for correct polarity description.

- ③ Replace the battery compartment cover.

BODY COMPOSITION ANALYSIS

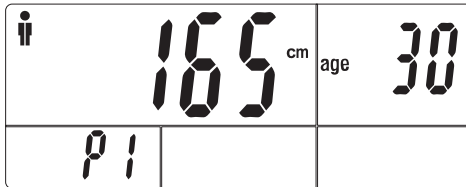
In addition to displaying the body weight, the device provides information on the percentage of bone, muscle, water and fat in the body based on the measurement and the entered user data (sex, age and height).

i The device is capable of storing data for 8 different users.

To enter personal data, proceed as follows:


SELECTING MEMORY STORAGE NUMBER

- 1 Press **SET** to enter the settings mode. The flashing **P1** icon indicates that the data of the first user can now be entered into the device memory.
- 2 Use **▲**, **▼** buttons to select storage location number.
- 3 Press **SET** to confirm your selection and enter the data of the next user.

SELECTING BODY TYPE AND GENDER

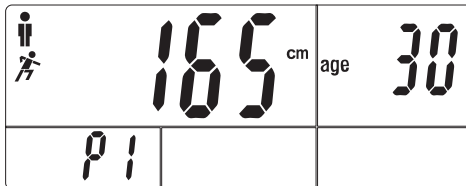- 1 Use arrows **▲**, **▼** to select your body type and gender represented by respective symbols. The following options are available: male **i**, male athlete **i/runner**, female **i** and female athlete **i/runner**.
- 2 Confirm your selection by pressing **SET**.

i To set the remaining parameters, proceed as above.

WARNING INDICATORS

Change battery – the display will show „Lo“. The appliance will switch off after a few seconds.

Overload – if the load exceeds 150 kg, the display will show „- - -“. Remove excess weight to avoid damaging the device.

! Do not exceed maximum load capacity of 150 kg.

Cleaning and maintenance

Clean with damp cloth. Do not submerge in water. Do not use abrasive detergents.

Ecology – Environment protection

Each user can contribute to protect the environment. It is neither difficult nor expensive. In order to do so:

Return cardboard packages to recycling points.

Throw polyethylene (PE) bags into appropriate containers.

Batteries contain substances harmful to the natural environment. Do not throw them away with other household waste. Return used batteries to recycling points. Before disposing of or returning the appliance to a recycling point make sure to remove the batteries from the battery compartment.

Return a worn out appliance to an appropriate collecting point as the components of the appliance may pose a threat to the environment.

Do not dispose into the domestic waste disposal!!!

Reference tables

BODY FAT PERCENTAGE

AGE	FEMALE				MALE			
	Underweight	Normal	Overweight	Obese	Underweight	Normal	Overweight	Obese
10~12	<12	12~22	22~30	>30	<8	8~18	18~24	>24
13~18	<15	15~25	25~33	>33	<8	8~18	18~24	>24
19~30	<20	20~29	29~36	>36	<8	8~18	18~24	>24
31~40	<22	22~31	31~38	>38	<11	11~20	20~26	>26
41~50	<24	24~33	33~40	>40	<13	13~22	22~28	>28
51~60	<26	26~35	35~42	>42	<15	15~24	24~30	>30
60+	<28	28~37	37~47	>47	<17	17~25	26~34	>34

BONE MASS

WEIGHT	FEMALE			MALE		
	<45 kg	45~60 kg	>60 kg	<60 kg	60~75 kg	>75 kg
BONES	1.8 kg	2.2 kg	2.5 kg	2.5 kg	2.9 kg	3.2 kg

BODY WATER

AGE	FEMALE	MALE
	Normal	Normal
10~12	>60	>64
13~18	>58.5	>63.5
19~30	>56	>62.5
31~40	>53	>61
41~50	>52	>60
51~60	>51	>59
60+	>50	>58

MUSCLE MASS




AGE	FEMALE	MALE
	Normal	Normal
10~99	>34	>40

The manufacturer/importer does not accept any liability for any damages resulting from unintended use or improper handling.

The manufacturer/importer reserves his rights to modify the product any time in order to adjust it to law regulations, norms, directives, or due to construction, trade, aesthetic or other reasons, without notifying in advance.

Notes

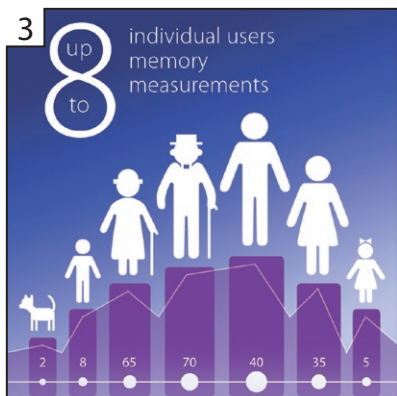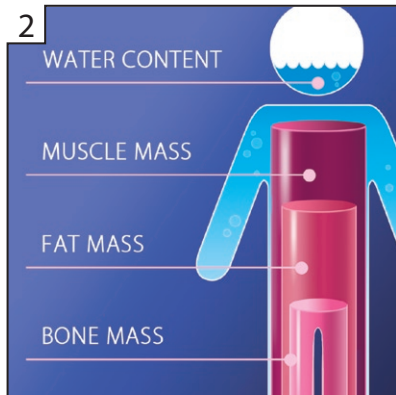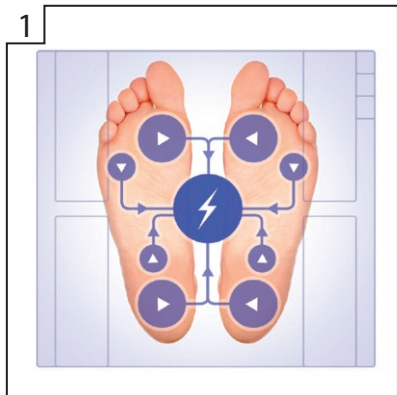A series of horizontal dotted lines for writing notes, consisting of 25 rows.

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EN 1. BIA (Bioelectrical Impedance Analysis) advanced through feet measurement system 2. Body composition analysis function – measures body bone, muscle, fat and water content 3. Measurements memory for up to 8 individual users (age, sex, height) 4. Maximum load – 150 kg (accuracy of measurement 100 g)